

HIV 1/2 Stat-Pak

For use with whole blood, serum, or plasma

Store Kits: 8 - 30° C

- Check kit before use. Use only items that have not expired or been damaged.
- Bring kit and previously stored specimens to room temperature prior to use.
- Always use universal safety precautions when handling specimens. Keep work areas clean and organized.


This outline is not intended to replace the product insert or your standard operating procedure (SOP).



1. Collect test items and other necessary lab supplies.



2. Remove device from package and label device with client identification number.

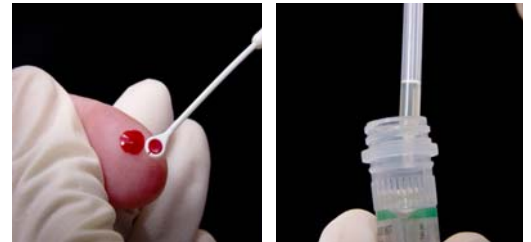

3. Collect approximately 5 µl of specimen using a new disposable loop or pipette.



4. Dispense the sample in the center of SAMPLE well.



5. Add 3 drops of buffer, holding vial vertically over the SAMPLE well.



6. Wait for 10 minutes before reading the results.



7. Read and record the results and other pertinent info on the worksheet.

HIV 1 / 2 Stat-Pak Test Results

Reactive

2 lines of any intensity appear in **both the control and test areas**.



Non-reactive

1 line appears in the **control** area and no line in the test area.



Invalid

No line appears in the **control** area. Do not report invalid results. Repeat test with a new test device even if a line appears in the test area.



